

Trinity™

Advanced Bearing Acetabular System



Acetabular system providing a seamless range of high performance bearings within a versatile shell featuring Biomimetic Cementless Technology.

Corin

Connected Orthopaedic Insight

Trinity™

Advanced Bearing Acetabular System


Patient matched solutions

Advanced implant design eases liner insertion and reduces risk of malalignment through an optimised taper¹, while a highly polished Prime™ rim aims to reduce soft tissue irritation.

Trinity-i™ shells are designed to maximise articulation size, for example matching a 36mm bearing to a 50mm shell. This means improved range of motion and joint stability can be achieved^{**}.

Shells are coated with vacuum plasma-sprayed pure titanium overlaid with a 20µm layer of calcium phosphate to maintain substrate roughness. This biomimetic cementless coating, which has over 20 years of clinical heritage^{2,3,4}, can promote long-term stability.

Trinity™ high performance bearings, including ECiMa™ (ultra-low wear⁵, vitamin-E blended polyethylene) and Dual Mobility, provide versatile solutions for a wide range of patients.

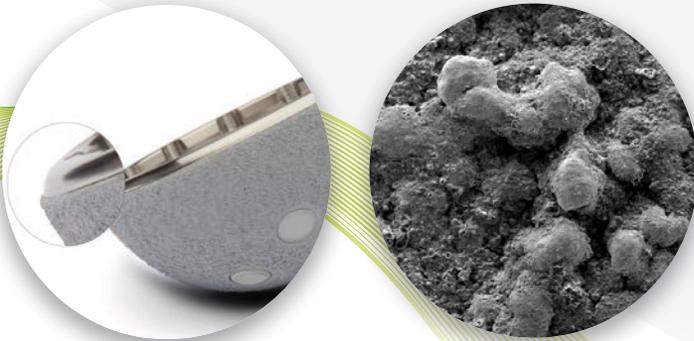

* Ceramic-on-ceramic bearing available in USA for evaluation only

**when comparing a 36mm bearing/50mm Trinity-i™ and a 32mm bearing and 50mm Trinity™.

1. Data on file at Corin Group Ltd.
2. Wood PLR, Deakin S. Total ankle replacement. The results in 200 ankles. J Bone Joint Surg (Br) 2003; 85-B:3:334.
3. 2. Saxler G, Temmen D, Bontemps G. Medium-term results of AMC unicompartmental knee arthroplasty. The Knee 2004; 11:39-355.
4. Schlueter-Brust KU, Kruse S, Bontemps G. Twelve year survivorship after cemented and uncemented medial unicompartmental knee arthroplasty. 15th EFFORT Congress June 2014.
5. Data on file at Corin Group Ltd.



coringroup.com

© 2019 Corin P No I1538 Rev01 12/2019 ECR2927